


AMRAD Medical



DR ELITE

IMAGE CAPTURE WORKSTATIONS

AMRAD® Medical DR Elite™ Image Capture Workstations

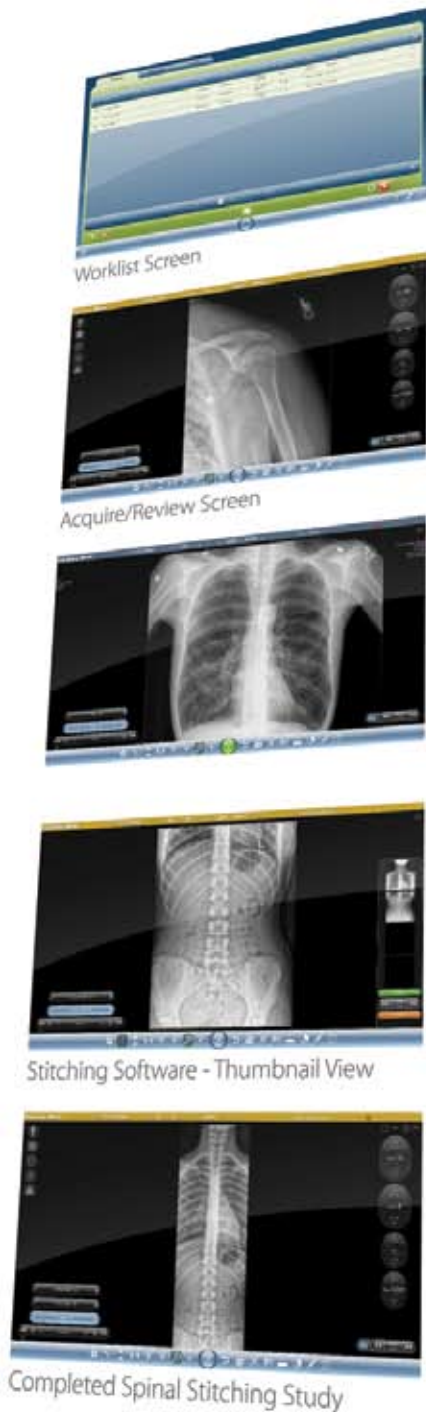

- Simple, clean, and intuitive two screen workflow: (Work list and Acquire/Review screen)
- Multi-touch or keyboard interface with touch gesture and animated sliding transparencies
- Image preview in 2 seconds from exposure
- Integration with Electronic Medical Record Systems (EMS) with full MWL, MPPS integrations with RIS/HIS providers
- Multiple automated exports that can be user configured with a single touch of a button
- Complete integration of dose area product meters (DAP) with all readings recorded for export as well as embedded in DICOM image tags

Advanced Features

- Generator Interface Control
- Auto-Positioning
- Low Dose Sensitivity Mode (LDSM) allows for as much as 50% dose reduction in dose sensitive studies while maintaining excellent diagnostic image results
- Rejection analysis with user customization of reports
- Embedded Operating System for network safe environments
- Embedded Help Functions for immediate support services
- Windows 10 Embedded Operating System

Options

- Stitching
- Workstations include 24/7/365 Support and 99% Uptime Guarantee



www.amradmedical.com

SUMMIT INDUSTRIES, LLC
7555 N. Caldwell Avenue, Niles, Illinois 60714
(p) 773-353-4000 | (p) 1-800-729-9729



LITR042-0820

SAI GLOBAL
ISO 13485
Quality