

ELEVATOR TABLE

The S211 elevator table features a 400-lb patient rating. It has foot-operated table controls (elevation, tabletop locks) and a standard 84x30 inch floating tabletop. A 76x30 inch tabletop is available.

While raising or lowering the tabletop it pauses at an adjustable working height.

The tabletop can only be mounted by sliding in on from the end.

TRAVEL AND POSITION SPECIFICATIONS

Tabletop travel	Longitudinal: ±15.5 in from center (84" Table Top) ±11 in from center (76" Table Top) Transverse: ±5 in
Tabletop height	Max 35 in Min 23 in
Patient load rating	400 lbs
Rear of lower table base to image center	Approximately 13.5 in

MOUNTING SPECIFICATIONS

Three mounting brackets are provided, each requiring a 3/8" floor anchor (*anchors not provided*). Select anchors appropriate for site floor construction. Recommended holding force per anchor: 500 lbs.


ELECTRICAL REQUIREMENTS

115 VAC @ 6.0A / 230 VAC @ 3.0A, 50/60 Hz; Standard 20-foot US line cord provided


SHIPPING INFORMATION

CONTENTS	DIMENSIONS	WEIGHT
Table base (crated)	57x33x39h in	550 lbs
Tabletop (crated)	88x35x4h in	80 lbs

Specifications subject to change without notice.

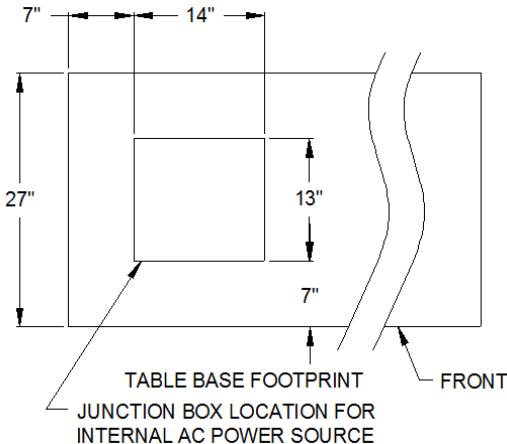

MOUNTING HOLE PATTERN

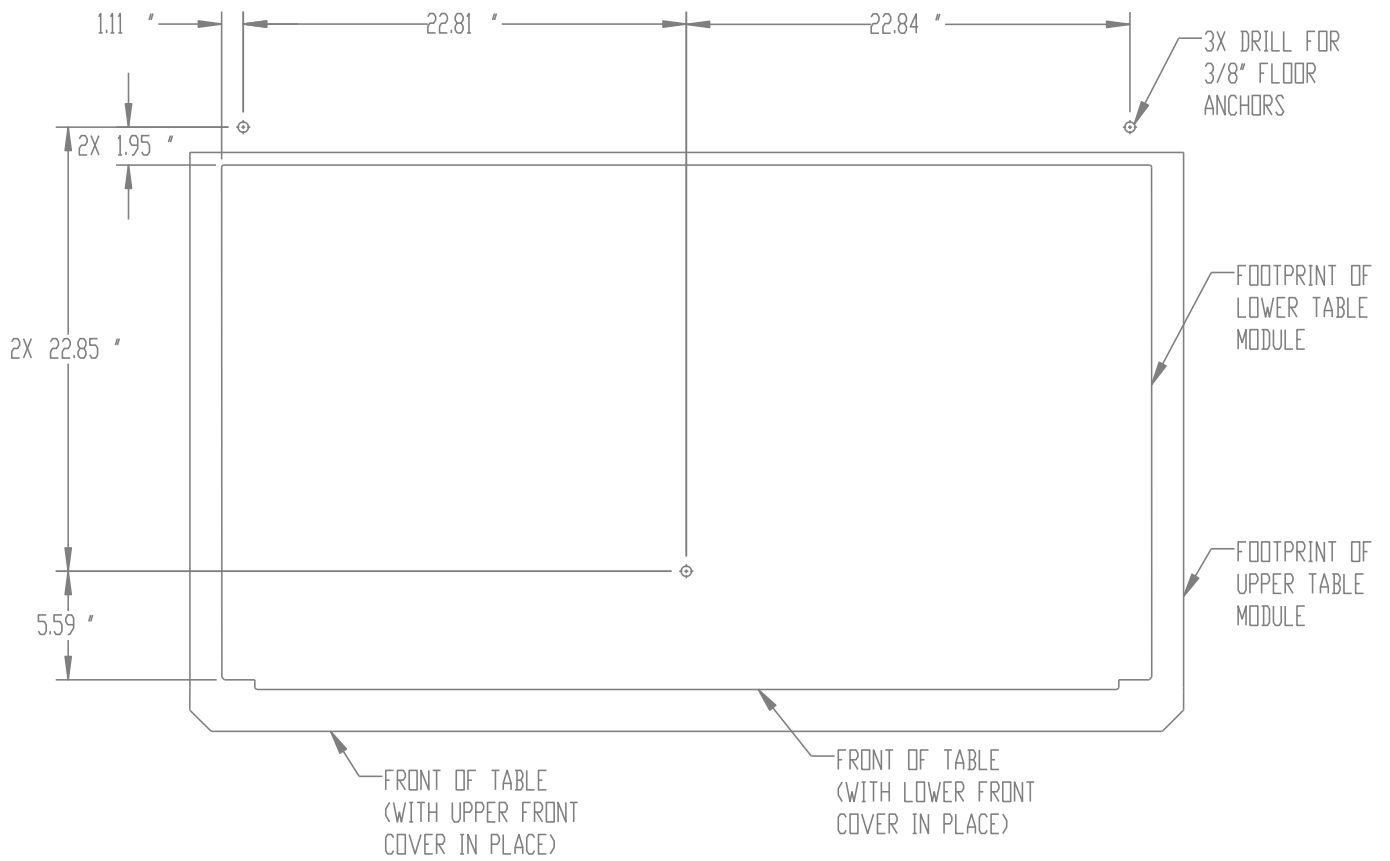

FOR BEST RESULTS, POSITION TABLE AND MOUNTING BRACKETS FIRST AND THEN USE MOUNTING HOLES IN BRACKETS TO DRILL FOR THE FLOOR ANCHORS.

ELEVATOR TABLE Model S211

