



S211 Elevating table

ELEVATOR TABLE

The S211 elevator table features a 400-lb patient rating. It has foot-operated table controls (elevation, tabletop locks) and a standard 84x30 inch floating tabletop. A 76x30 inch tabletop is available.

While raising or lowering the tabletop it pauses at an adjustable working height.

The tabletop can only be mounted by sliding in on from the end.

TRAVEL AND POSITION SPECIFICATIONS

Tabletop travel

Longitudinal: ± 15.5 in from center

Transverse: ± 5 in

Tabletop height

Max 35 in

Min 23 in

Patient load rating 400 lbs

Rear of lower table base to image center Approximately 13.5 in

MOUNTING SPECIFICATIONS

Three mounting brackets and mounting locations provided, each requiring a 3/8" floor anchor (*anchors not provided*). Select anchors appropriate for site floor construction. Recommended holding force per anchor: 500 lbs.

ELECTRICAL REQUIREMENTS

115 VAC @ 6.0A / 230 VAC @ 3.0A, 50/60 Hz; Standard 20-foot US line cord provided


SHIPPING INFORMATION


CONTENTS DIMENSIONS WEIGHT

Table base (crated) 57x33x39h in 550 lbs

Tabletop (crated) 88x35x4h in 80 lbs

Specifications subject to change without notice.

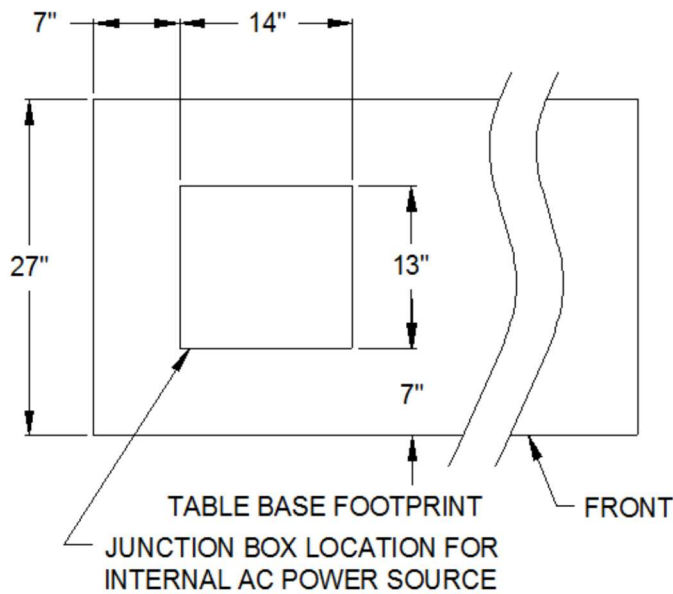

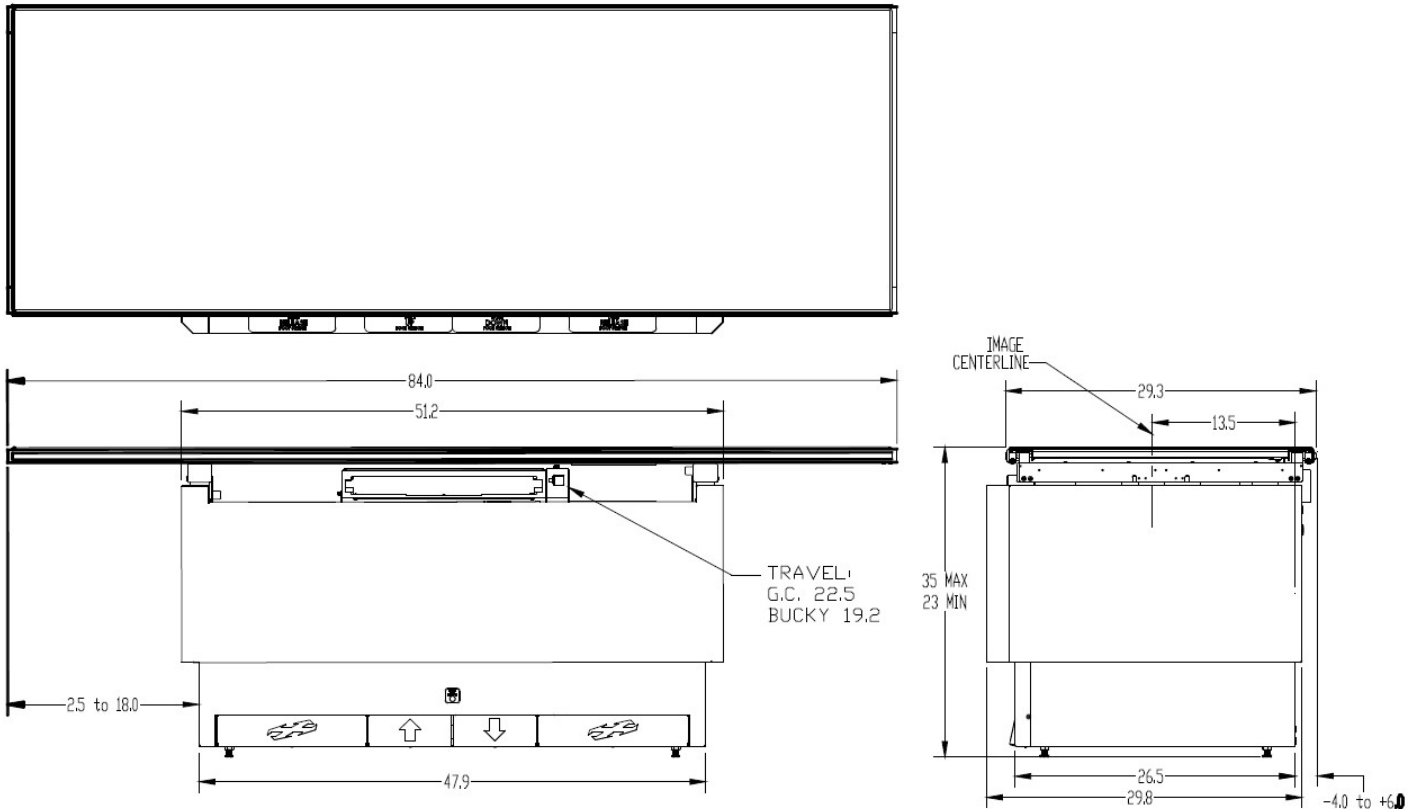

AMRAD Medical

SUMMIT INDUSTRIES, LLC
7555 N. Caldwell Avenue, Nilis, Illinois 60714
(p) 773-353-4000 | (p) 1-800-729-9729