



FWFC Classic

FLOOR-TO-WALL / CEILING TUBE STAND SYSTEM



The **FWFC Classic System** can be configured to meet the needs of any department. This full-featured floor-to-wall / ceiling system is available with either 8 ft. or optional 10 ft. rail to perform weight-bearing studies. Plus the FWFC Classic System comes with an industry-leading 5-year parts warranty.

FWFC Classic

Floor-to-Wall /Ceiling Tubestand

- Cross-table exams are easily performed using the +/- 180° rotational and +/- 5 in. transverse positioning capabilities
- Floor-to-wall or floor-to-ceiling mount design easily fits into 8' ceiling height (93" ceiling minimum)
- Available with either 8 ft. or optional 10 ft. rail
- Certified 150 kVp manual collimator with swivel mount
- Laser positioning light
- LED light field and tape measure



A compact design for facilities with space constraints that is both economical and versatile.



DR UPGRADEABLE
ASK ABOUT OUR DR UPGRADE SOLUTIONS



Four-way Float Top Table

- Four-way, non-elevating floating top permits easy transfer and positioning of patients
- Recessed foot-treadle allows for safe patient transfer
- 28 inch longitudinal travel; 10 inch transverse travel



Positioning features such as the comfortable patient chin rest and one-touch release handle allow the technologist to focus on the patient.



Wallstand

- Designed for ease of use and simplicity, the wallstand lowers for weight bearing exams like standing knees
- Ergonomically designed handle allows for easy positioning of the wall receptor



Chin rest



One-touch release handle

Specifications

FWFC Classic Tube Stand (S225)

- Floor-to-wall or floor-to-ceiling mount with 8 ft. rails
- Rotational base, +/- 360° with 90° detents
- Transverse tubearm +/- 5 in., platform tube mount
- 14 in. - 76 in. vertical travel
- Angulation dial, operator handgrips, and electric locks

FWFC Classic Wall Stand (S109)

- Counterbalanced column, floor-to-wall mounted
- Grid Cabinet, 17 in. x 17 in.
- Grid: 10:1 ratio, 103 LPI
- Heavy-duty deluxe cassette tray 06338
- Specify left or right load
- Electric locks
- Accepts all size detectors

Non-Elevating Table (S210)

- Four-way float top
- 84 in. x 30 in. table top, table height 31.5 in. (Optional 76-inch top is also available)
- 400 lb. rated patient load
- Heavy-duty deluxe cassette tray


Generator options

- 30 - 50kW generator options

Specifications are subject to change without notice



is a Trusted Brand of Summit Industries

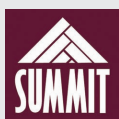

Since 1984

SUMMIT INDUSTRIES, LLC
7555 N. Caldwell Avenue, Niles, Illinois 60714
(p) 773-353-4000 | (p) 1-800-729-9729
www.summitindustries.net



LITX005-0323
FWFC Classic