


Agility



AGILITY

Mobile C-Arm
All-in-One Fluoroscopy System

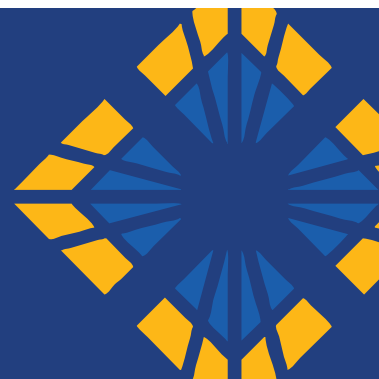

Simple Unique

Usability Comes First

DISPLAY UNIT

Balanced articulated arm. Reduced dimensions, yet still large image display with its 27" landscape monitor. Touchscreen interface for adapting and post processing images. Wide sideways extension of the articulated arm with height and angle adjustable monitor.

Three-way articulated monitor support arm for best frontal vision. The monitor position can be adjusted from both sides in height, angle and tilt.



FLAT PANEL

Available acquisition surfaces range from 21x21cm to 30x30cm for high resolution images. Removable anti-scatter grid for low dose applications.



MODULATED POWER

X-ray emission in pulsed mode with modulated power to suit every procedure type. Excellent dose/image quality ratio.





REMOTE DISPLAY

High resolution WIFI video transmission for remote display of the images on either an optional wireless monitor cart, monitor ceiling suspension or wall-mount monitors.



ABSOLUTE OPERABILITY

Swiveling multi-touch console allows full access from any side. The operator benefits continuously from an unobstructed view of the live images, while having the controls always within his reach.

Optional integrated printer for final procedure report output.



COMPACT DESIGN

Large diameter, lightweight mobile C-arm for easy access to the operating table. Operators will appreciate the effortless, fast and accurate positioning around the patient. The single unit, compact design makes it easy to move within and between operating rooms.



Ready to Go

Perfectly balanced orbital and angular rotation of the C-arm for fast and precise positioning. Responds perfectly to operator touch and feel.

DISPLAY

View images with a maximum comfort. The monitor height and tilt are adjustable thus keeping it always in front of the operator. No obstacles on the floor.



UP/DOWN MOVEMENT

The large vertical movement range of 18 in (450 mm) turns this C-arm into a very versatile unit.




ORBITAL ROTATION

Vast 160° orbital rotation range 160° (from +65° to -95°) speeds up angular and lateral projections without having to move the stand around the operating table.



ANGULAR ROTATION

Precise movements and wide angles of more than +/- 210° around and in front of the patient.

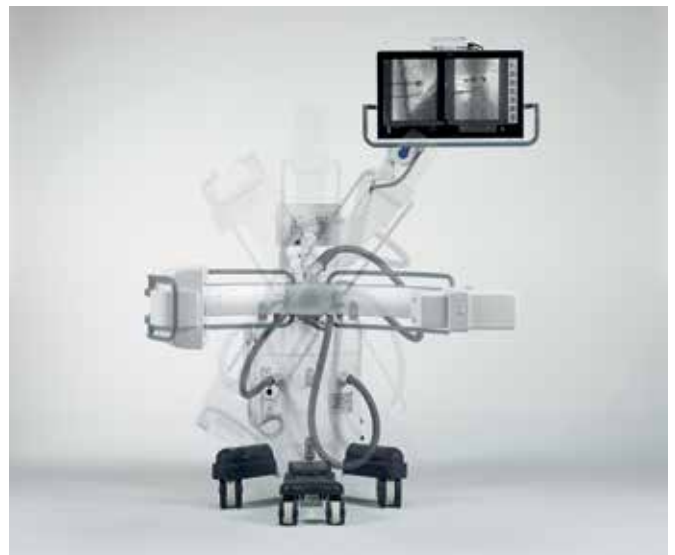

360° Vision

DISPLAY UNIT

The ultra-wide 27" monitor allows for simultaneous display of live and stored images. Thanks to the 4 MP resolution, the monitor shows even the smallest details. The wide travel of the articulated arm lets you move the monitor comfortably into the ideal position for use or visualization.



ADDITIONAL DISPLAY UNIT

The optional wireless monitor cart features WIFI video transmission, powered by a strong lithium-ion battery for increased autonomy. Wide field of view, duplicating the main monitor.



CONSOLE

High resolution, high brightness console with a complete operation panel providing all tools for acceptance, selection, edition, movement, transmission, compression, reduction, enlargement... a simple touch is enough.

Quality Features

Removable grid for low dose pediatric surgery and gridless procedures.



X-ray emission controlled by either a foot switch or a multi-function hand switch.



Display certified for medical use and in conformity with the relevant standards for diagnostic imaging.



The best possible connectivity: USB, Ethernet, HDMI, and WiFi.





Zero dose centering thanks to laser localizers on X-ray tube and on flat panel detector side.



Optional wireless footswitch for more freedom and less cables on the floor.



Instant login via the optional NFC device.



STILLE Medstone Tables

Hospital grade surgical imaging tables by STILLE offer superior versatility. Several models are available for purchase.

Model dependent features include

- Tables with up to 5 motorized movements, including:
 - Height Adjustment
 - Trendelenburg
 - Lateral tilt
 - Longitudinal travel
 - Lateral travel
 - Plush 2" thick mattress pad with Sofcare™ technology
 - 500 lb. patient weight capacity
 - Infection control covers for several models.
- See individual model specifications for details.



PM TABLES

Features a contoured table top with a facial cutout. Ideal for orthopedics, pain management and electrophysiology.

Table width: 22 in



ERCP TABLES

Ideal for any procedure that requires a wide table top and heavier load capacity.

Table width: 26 in



CV TABLES

Free float technology allows for head-to-toe imaging with accuracy, speed and freedom of movements to improve efficiency in the OR. Ideal for procedures such as cardiovascular, endovascular and fluoroscopic.

Table width: 22 or 24 in (dependent on model)

Contact your Summit sales representative to determine which table will best fit your needs.

Support

Anytime, Anywhere



GLOBAL SUPPORT

Summit Industries boasts a certified quality management system since 1994 and was one of the first medical device manufacturers in Europe to do so. Summit Industries provides technical support for its products throughout their useful life, wherever they may be used.



REMOTE SUPPORT

Connectivity is native to the AGILITY, making remote support as easy as can be, whether you need to perform efficiency tests, software upgrades or maintenance.



LOCAL SUPPORT

Summit Industries runs its own training academy for personnel with a mission to provide technical and application support to Summit Industries products.

Technicians from all over the world are trained by Summit Industries, receiving a certificate of qualification to enable them to service AGILITY.

PIONEERS IN RADIOLOGY SINCE 1984

A NETWORK OF SERVICE + SUPPORT

Enjoy all the advantages of our successful service network. From the first shot to the last managed image, we work to continually improve your imaging performance and patient experience.

Summit Industries' successful service network surrounds you with a dynamic team of exceptional experts. Your organization and your patients will benefit from the experience and best practices that only Summit Industries and its authorized dealer network can offer.



AMRAD Medical | AGILITY

SUMMIT INDUSTRIES, LLC
7555 N. Caldwell Avenue, Niles, Illinois 60714
(p) 773-353-4000 | (p) 1-800-729-9729

amradmedical.com



LITR053-0623