

# *SIMPLICITY*



## **SIMPLICITY** FW / FC SYSTEM

FLOOR-TO-WALL / FLOOR-TO-CEILING  
MOUNTED TUBE STAND SYSTEM

# Compact and affordable

The AMRAD Medical **SIMPLICITY** FW/FC System is a compact configuration designed for medical imaging centers and urgent care facilities looking to maximize space and minimize room modifications. It features a floor rail and an upper rail that can mount to the wall or ceiling to accommodate longitudinal travel for quick and accurate positioning.

The simple, yet rugged design makes this an affordable radiographic system built for longevity and performance making it an ideal choice for lower volume facilities.



## PERFORM A FULL RANGE OF EXAMS!

- ✓ Table Exams
- ✓ Upright Chest
- ✓ Extremities
- ✓ Weight Bearing Studies

## INCREASE EFFICIENCY!

The AMRAD Medical **SIMPLICITY** is designed with our Easy-Read Angulation Dial for quick and accurate identification and adjustment of angles.

Specially designed detents of 40" to table and 40" and 72" to wall make positioning quick and efficient.

## RADIOGRAPHIC TABLE

The four-way float top table allows for comfortable and easy patient positioning with its 31" access height. The recessed foot-treadle allows for hands-free table operation as well as patient safety. The table has a 400 lb. rated patient load.



## WALL STAND

Designed for flexibility, the **SIMPLICITY'S** wall stand has a maximum vertical travel range of 72<sup>3</sup>/<sub>4</sub>" and can go as low as 9<sup>3</sup>/<sub>4</sub>" for weight bearing exams like standing knees.

The One-Touch Release Handle allows for easy positioning of the wall receptor.



*Comfortable chin rest*



*Push button Lock Release for quick positioning*

## DR Ready



All AMRAD Medical systems are DR compatible. Ask about our DR Packages available with 100 micron Amorphous Silicon (a-Si) TFT CsI (Cesium Iodide) detectors that ensure superior image quality at a lower dose.

# Specifications

## Radiographic Table (S210)

### FOUR-WAY FLOAT TOP TABLE

- Four-way float top table with foot treadle control
- 400 lb. rated patient load
- +/- 5 in. transverse travel, 28 in. longitudinal travel
- 84 in. x 30 in. table top, table height 31.5 in.
- 17 in. grid cabinet
- Grid: 10:1 ratio, 103 LPI
- Heavy-duty deluxe cassette tray

## Wall Stand (S109)


### WALL-MOUNTED WALL STAND (17"X 17")

- Fully counterbalanced
- 17 in. grid cabinet with front panel travels from floor to 76 in.
- Grid: 10:1 ratio, 103 LPI, specify focal range
- Specify left or right load
- Heavy-duty deluxe cassette tray 06338
- Electric locks
- Operates on 24VDC

## Tube Stand (S225)

### FLOOR-TO-WALL/FLOOR-TO-CEILING MOUNTED TUBE STAND

- Transverse tubearm +/- 5 in., platform tube mount
- Rotational base, +/- 360° with 90° detents
- Fixed length angulating tube arm
- Easy to read angulation dial
- Vertical travel: 13.5" minimum to 76.8" maximum (focal spot to floor)
- Electric locks: vertical, horizontal, and tube arm angulation, all locks release
- Fully counterbalanced
- 8' or optional 10' floor rail
- Operates on 24VDC powered by Nexus box

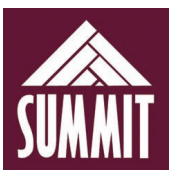

### A NETWORK OF SERVICE + SUPPORT

*Summit Industries' successful service network surrounds you with a dynamic team of exceptional experts. Your organization and your patients will benefit from the experience and best practices that only Summit Industries and their Authorized Dealer Network can offer.*

PIONEERS IN RADIOLOGY SINCE 1984

AMRAD Medical | SIMPLICITY

### SUMMIT INDUSTRIES, LLC

7555 N. Caldwell Avenue, Niles, Illinois 60714  
(p) 773-353-4000 | (p) 1-800-729-9729



NSF-ISR

Registered to ISO 13485



LITR033-0124