




OTS ELITE

CEILING MOUNTED TUBE SUPPORT
AUTO TRACKING SYSTEM



Performance, Versatility, and Low Cost of Ownership

The AMRAD Medical OTS Elite delivers maximum capability and is designed for use in high-volume orthopedics, imaging centers, free-standing ERs and hospitals.

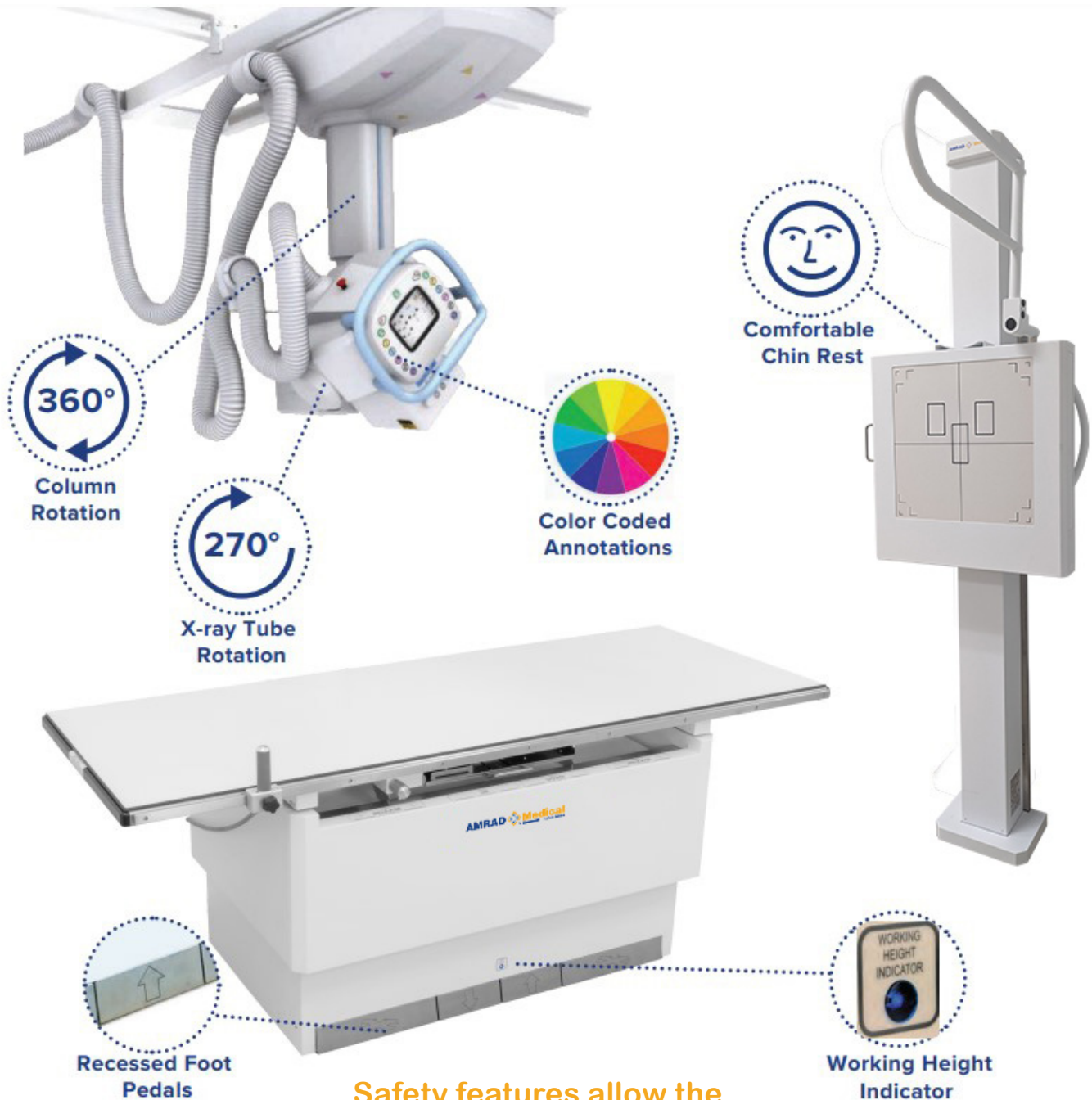

Safety features allow the technologist to focus on the patient





M201 AUTO-TRACKING OVERHEAD TUBE

Ceiling-mounted tube with auto-tracking to the table and wall stand that can be turned on or off at the touch of a button.

- ✓ One-touch lock release and color-coded annotations make positioning fast, improving patient throughput
- ✓ Vertical axis is motor driven for ease of movement
- ✓ Electromagnetic locks with safety ON in case of power outage

S222 RADIOGRAPHIC TABLE

The elevating, four-way float, flat top allows for easy transfer and positioning of patients.

- ✓ 650 pound patient capacity
- ✓ Recessed four-pedal design provides easy operation and patient safety
- ✓ Collision protection and fail-safe locks for added patient safety



Table Elevation Control Handle (T.E.C.H.) with lock release

J1000 WALL STAND

The fully counterbalanced wallstand with rigid construction ensures precise alignment.

- ✓ Center or side-mounted for patient and/or configuration flexibility
- ✓ Side-mounted and overhead hand grips for patient stability and positioning
- ✓ Collimator light can be activated from the wallstand, improving productivity



Ergonomic lock release handle on wallstand for ease of positioning

OTS Elite Specifications

Radiographic Table (S222)

ELEVATING FOUR-WAY FLOAT TOP TABLE

- Foot treadle operation of table top and height movements
- 650 lb. rated patient load with collision protection
- Variable table height, 23 in. to 34 in.
- +/-5 in. transverse travel, 30 in. longitudinal travel
- 84 in. x 30 in. flat, seamless table top
- Optional 35.5 in. width table top is also available

Wall Stand (J1000)

WALL-MOUNTED WALL STAND (17"X 17")

- Fully counterbalanced
- 17 in. x 17 in. grid cabinet with front panel travels from floor to 76 in.
- Grid: 10:1 ratio, 103 LPI, specify focal range
- Specify left or right load
- Heavy-duty deluxe cassette tray 06338
- Electric locks
- Operates on 24VDC
- Optional Long 17 x 51 in. format wall stand is also available

Overhead Tube Support (M201)

CEILING MOUNTED TUBE SUPPORT

Recommended ceiling height 9' to eliminate any positioning challenges

- 15 feet / 4.6 meter longitudinal rails – longitudinal travel 192 in.
- 10 feet / 3 meter bridge – transverse travel 95 in.
- Motorized vertical column – vertical travel 63 in.
- Tube column rotation - 180° with detents every 45°
- Tube angulation $\pm 135^\circ$ with detents every

Generators

- Choose from our S Series, E Series or M Series 40 kW - 80 kW HF Generator options



Both our OTS Elite and OTS Classic systems offer equal features and earn equal reimbursement as that of the more costly competition . . . all of this at 5-times the warranty!



SUMMIT INDUSTRIES
A NETWORK OF SERVICE & SUPPORT

Summit Industries' successful service network surrounds you with a dynamic team of exceptional experts. Your organization and your patients will benefit from the experience and best practices that only Summit Industries and their Authorized Dealer Network can offer.

AMRAD  **Medical**
by **Summit Industries**

SUMMIT INDUSTRIES, LLC

7555 N. Caldwell Avenue, Niles, Illinois 60714

(p) 773-353-4000 | (p) 1-800-729-9729



LITR043-1124

www.amradmedical.com