

Radiographic Floating Top Tables



A4000 Elevator Table

AMERICOMP

All the Features You Need at a Price You Can Afford


Since 1987, the quality standards of Americomp have been unsurpassed. Americomp radiographic components remain the affordable choice in today's marketplace, offering the same long-term dependable performance that has become synonymous with the Americomp name. Our manufacturer's parts warranty secures your investment today and guarantees your investment value tomorrow.

Float Top Tables

The A4000 elevating float top table provides an elevation range of 22 1/2" to 34" (57 cm to 86 cm) for effortless patient transfer. It is very well suited for immobile, pediatric or geriatric patients. The smooth, quiet lift mechanism and high-capacity motor safely handles patients weighing up to 400 lbs. (182 kg).

The A5480 float top table is 31" (79 cm) in height with a welded steel pedestal base that facilitates patients weighing up to 400 lbs. (182 kg).

Both tables feature a generous table top workspace, allowing for comfortable and easy patient positioning. Precision stainless steel ball bearings float the table top in four directions. The heavy-duty table top locks provide a stable surface for fast and safe transfer of patients even during power outages. The recessed foot controls within the table base permit hands-free operation, while restricting inadvertent patient release.

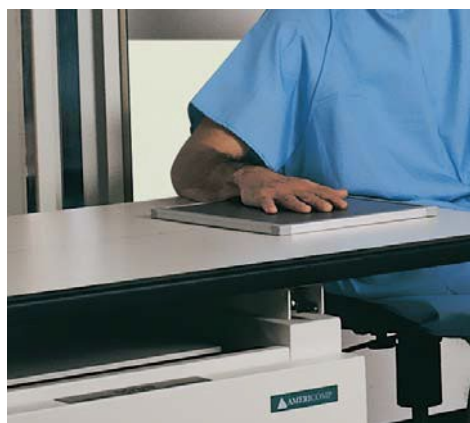

The large, flat working surface of the table provides images without distortion.



A5480 Four-way Table

High Performance Generators

Complement any Americomp radiographic table with a high-quality, state-of-the-art generator. Contact your Americomp dealer for details.

Digital-ready

Americomp radiographic components are completely compatible with most digital imaging systems. Digital imaging technology may be added at any time.

Americomp

7555 N. Caldwell Ave.
Niles, Illinois 60714
1-800-458-1782
www.americompus.net